





Electric Bed

for a safer healing environment



Feature

Angle Indicator

Medical personnel can gather clear and easy readings of the back raise angle of the patient. The angle indicator assists medical personnel to raise the backs of patients above 30° to prevent ventilator-associated pneumonia (VAP).



SSR (Stress-Saver Raise)™

Stress on the patient's body is reduced during back raise, as the back section of the bed frame extends to accommodate a natural sitting up process. This also prevents the patient from sliding down towards the foot of the bed.





Fall Prevention

Safety Side-rails

Practical features such as side-rails that move along with the back section of the bed minimize the risk of patients falling when resting in a flat or up right position.

(complies with IEC60601-2-52)



Low Bed Height 350mm*

A low height bed 350mm* ensures that patients can firmly place their heels on the floor while the ergonomic side-rail design provides patients with additional support when standing up.


 * Height from floor to mattress base surface



Easy to Use

Compact Design

Other than providing a comfortable space for patients to rest on, the compact design of the overall bed width facilitates greater ease when maneuvering the bed around sharp corners, narrow entrances and elevators. As part of its compact design, gaps between the bed and its components are reduced. This allows medical personnel to closely position a stretcher next to the bed for quick and easy transfer of patients.

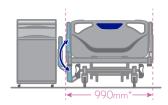

*When the side-rails are lowered: 970mm

Specifications

	MA-54380
Total Length	2,200mm
Total Width	990mm (when the side-rails are lowered: 970mm)
Total Height	755 to 1,105mm
Mattress Base Height*	350 to 700mm
Mattress Size	2,000 x 860mm
Safe Working Load**	230kgf (2,250N)
Caster	125mm dia. single-wheel casters with total locking system
Back Raise	0° to 70°
Knee Raise	0° to 25°
Trendelenburg/ Reverse- Trendelenburg	0° to 12°

PARAMOUNT BED VIET NAM TRADING CO.,LTD

HMC Building, 193 Dinh Tien Hoang street, DaKao Ward, Dist 1, Hochinminh city, Vietnam

Tel: +84 28 38 206 078 Fax: +84 28 38 206 079

https://www.paramount.co.jp/english/



- $^{\star}\,$ Indicates the height from the floor to the mattress base surface.
- ** Maximum load with which the bed can be operated.

 (The total load of user's weight and incidental items including the mattress and options.)

 $Battery, X-ray\ imaging\ base\ and\ Individual\ locking\ casters\ are\ also\ available.$ For more information, please ask Paramount\ Bed\ or\ your\ distributor.